



The Signature System and the Vanguard Total Knee

- Offers 10 femoral and 9 tibial sizes with complete interchangeability providing custom fit for the largest percentage of the patient population
- Designed to allow 145° of flexion with no additional bone resection
- Narrow anterior cortex to avoid overhang in smaller sized femurs
- Patella friendly design with a long deep trochlear groove, without the need of additional bone resection
- Thin, rounded femoral edges to allow pain free articulation with the soft tissues
- 3 levels of Cruciate retaining and 2 levels of Posterior stabilised inserts

This publication has been issued by:

European Central Marketing
Waterton Industrial Estate
Bridgend, South Wales
CF31 3XA, United Kingdom
Tel: +44 (0)1656 655221
Fax: +44 (0)1656 645454

Authorised Representative in the EU

Biomet UK Ltd.
Waterton Industrial Estate
Bridgend, South Wales
CF31 3XA, United Kingdom
Tel: +44 (0)1656 655221
Fax: +44 (0)1656 645454

Disclaimer

This publication and all content, artwork, photographs, names, logos and marks contained in it are protected by copyright, trademarks and other intellectual property rights owned by or licensed to Biomet or its affiliates. This brochure must not be used, copied or reproduced in whole or in part for any purposes other than marketing by Biomet or its authorised representatives. The use for any other purposes is prohibited.

Biomet does not practice medicine and is not responsible nor does it recommend for the selection of any particular orthopaedic implant or surgical technique for use on a specific patient. The surgeon who performs any implant procedure is responsible for determining and utilising the appropriate techniques for implanting prosthesis in each individual patient.

© Biomet 2011

Signature is a registered trademark from Biomet UK Ltd.



One Surgeon. One Patient.



Signature

Personalised Arthritis Care

BIOMET[®]


Signature Personalised Arthritis Care

Surgeon designed. Patient matched.

The Signature System and the Vanguard Total Knee

Combining the precision of patient specific positioning guides developed from MRI with the most flexible knee on the market allows for a custom-fit of the femoral and tibial components independently to address the largest possible percentage of the population.

- All of the alignment and accuracy benefits of navigation without the additional OR time, fixation of the arrays or capital investment
- MRI Scan leads to patient specific plan of implant size and positioning
- Surgeon has ability to alter the plan before approval of the Signature guides and their manufacture
- In the OR the patient specific Signature guides find their unique position on the patients anatomy, providing exact pin placement for the distal femoral resection guide, 4-in-1 cutting block and the proximal tibial cutting guide

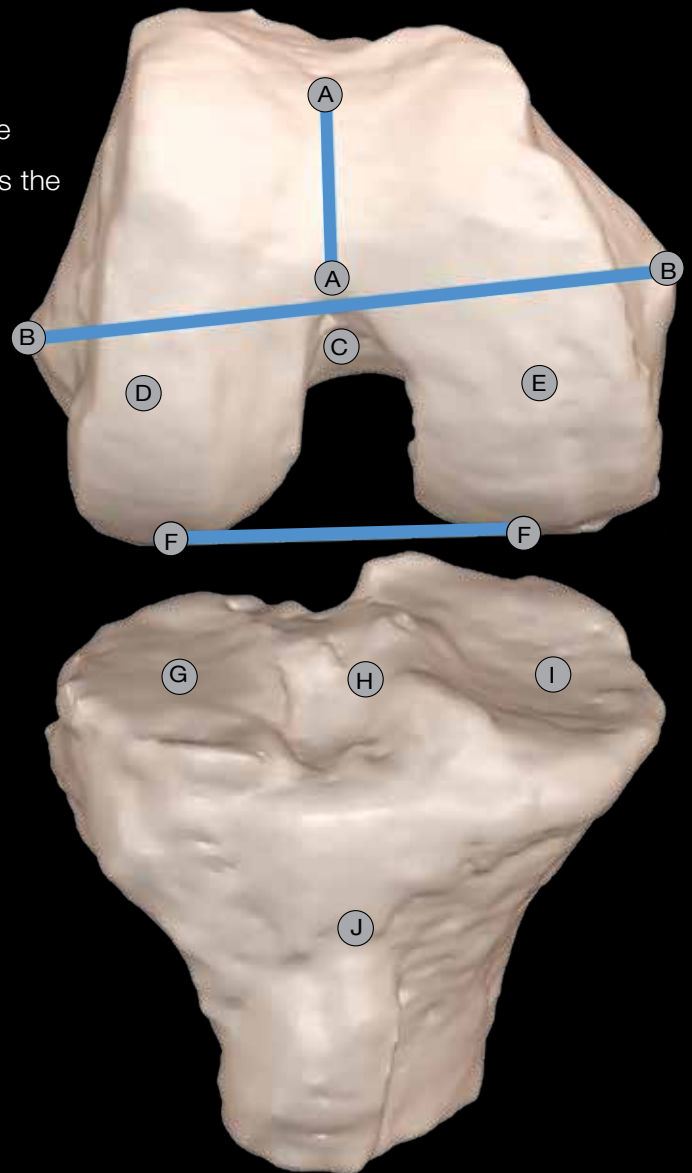

Signature Planning Landmarks

- A Anterior/Posterior Axis
- B Epicondylar Axis
- C Distal Femoral Mechanical Axis Landmark
- D Lateral Distal Femoral Condyle
- E Medial Distal Femoral Condyle
- F Posterior Condylar Axis
- G Lateral Plateau
- H Proximal Tibial Mechanical Axis Landmark
- I Medial Plateau
- J Medial One-third of the Tibial Tubercle

Signature

Workflow Timeline

1. MRI
Validation

2. MRI
Scan

Signature User Interface

- Data from the patient MRI generates a preoperative plan shown below
- The user friendly planning software allows the surgeon to fine tune the preliminary plan
- Final guides will incorporate the surgeons saved preoperative alterations



3. Upload
Image Data

4. Schedule
Surgery

5. Download
Initial Plan

6. Implant
Planning

7. P
M

Signature Positioning Guides

Femoral Positioning Guides Establish:

- Distal femoral resection level (9mm default)
- Distal femoral resection
 - Valgus angle (perpendicular to femoral mechanical axis default)
 - Flexion angle (3 degree default)
- Femoral component rotation (epicondylar axis default)
- Anterior/posterior position (9mm posterior condylar resection default)
- Femoral component sizing



Tibial Positioning Guides Establish:

- Proximal tibial resection level (10mm default)
- Proximal tibial resection
 - Varus/valgus angle (perpendicular to tibial mechanical axis default)
 - Tibial slope angle (3 degree default)